

KT I Lab Course

Angular Correlation

Assistant:

Alexander Kish

Students:

Fabian Stäger, Manuel Sommerhalder

Date of experiment:

24.01.2018







Contents

1	Introduction	3			
	1.1 The ⁶⁰ Co Decay	3			
	1.2 Angular Correlation	3			
	1.3 Experimental Concept	4			
2	Measurement Principle	5			
	2.1 Experimental Set-Up	5			
	2.2 Electronic Modules	7			
	2.3 Scintillation Counters and Photomultipliers	9			
3	Energy Spectrum	10			
4	Measurement Log 12				
5	Data Analysis	14			
	5.1 Time Calibration	14			
	5.2 Gaussian Fit	15			
	5.3 Angular Correlation Function	16			
6	Conclusion	18			
	6.1 Results	18			
	6.2 Discussion of Error Sources	18			
7	Appendix	20			
	7.1 Energy Spectrum	20			
	7.2 Measurement Log	21			
	7.3 Data Analysis	22			
	7.4 References	22			



1 Introduction

The goal of this experiment is to measure the angular correlation function of γ -rays in a 60 Co decay and compare it to theoretical predictions.

1.1 The ⁶⁰Co Decay

⁶⁰Co is a synthetic isotope of cobalt with a half-life of 5.1 years. Through β^- decay it disintegrates into an excited ⁶⁰Ni atom with an energy of 2505 keV and l=4 with positive parity. This nickel isotope thus emits two successive gamma rays with energies of 1172 keV and 1332 keV in order to reach its stable state of 0⁺. Since the intermediate state (2⁺) has a lifetime of around 1 ps, the two γ-rays can due to the finite experimental time resolution be treated as coincident. A schematic diagram of the ⁶⁰Co decay is provided in figure 1.

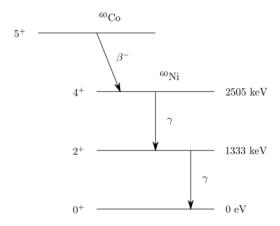


Figure 1: Schematic illustration of the ⁶⁰Co Decay. [1]

1.2 Angular Correlation

The two successive γ -rays are predicted to be emitted anisotropically relative to each other according to the general angular correlation function [1]:

$$W(\theta) = 1 + \sum_{i=1}^{l} a_i \cos^{2i} \theta \tag{1.1}$$

This function describes the probability for the latter gamma ray to be emitted at an angle θ from the first one. 2l is the order of the lowest multipole in the cascade. In



case of the ⁶⁰Co decay, both gamma rays are emitted during a transision of $\Delta l = 2$ and positive parity, so the dominant contribution is the electric quadrupole and the correlation function reduces to:

$$W(\theta) = 1 + a_1 \cos^2 \theta + a_2 \cos^4 \theta \tag{1.2}$$

The coefficients a_1 and a_2 have been theoretically predicted by Dr. Hamilton [2] for all combinations of angular momenta involved by consideration of state transitions using the Clebsch-Gordan-coefficients. These are summarized in figure 2

J ₂	J_2	Multipoles	a_1	az
1	0	Dipole-Dipole	1	0
1	1	Dipole-Dipole	-1/3	0
1	2	Dipole-Dipole	-1/3	0
1	1	Quadrupole-Dipole	-1/3	0
1	2	Quadrupole-Dipole	3/7	0
1	3	Quadrupole-Dipole	-3/29	0
2	3	Dipole-Quadrupole	-3/29	0
2	2	Dipole-Quadrupole	3/7	0
2	1	Dipole-Quadrupole	-1/3	0
2	0	Quadrupole-Quadrupole	-3	4
2	1	Quadrupole-Quadrupole	5	-16/3
2	2	Quadrupole-Quadrupole	-15/13	16/1
2	3	Quadrupole-Quadrupole	0	-1/3
2	4	Quadrupole-Quadrupole	1/8	1/24

Figure 2: Coefficients a_1 and a_2 for different combinations of J_1 and J_2 for an atom deexciting into the ground state $J_1 = 0$ [1]

Since the angular momenta of $J_3 = 4$, $J_2 = 2$, $J_1 = 0$ are assumed for the ⁶⁰Co decay, the expected correlation function has the following explicit form:

$$W(\theta) = 1 + \frac{1}{8}\cos^2\theta + \frac{1}{24}\cos^4\theta \tag{1.3}$$

1.3 Experimental Concept

The experimental setup, more throroughly described in section 2.1, consists of two scintillation counters that detect the coincident γ -rays from different and adjustable relative angles. The rate of detected events is then plotted and compared to the theoretical description from formula 1.3.



Measurement Principle 2

Experimental Set-Up 2.1

The 60 Co-sample is placed in the center of a goniometer table where a large scintillation counter D_l is set up at a fixed position, measured to be at a distance of $r_l = 216$ mm from the sample, pointing to the center. A smaller scintillation counter D_s is mounted on a rotatable arm, also pointing to the center. The angle relative to the fixed detector can be read off the edge of the goniometer table. The distance r_S from the small detector to the sample can be adjusted. Here we chose 87 mm in order to achieve the same solid angle coverage as the large detector. This was deduced from the definition of the solid angle:

$$\Omega = \frac{A}{r^2} = \frac{\pi}{4} \frac{d^2}{r^2} \tag{2.1}$$

$$\Omega = \frac{A}{r^2} = \frac{\pi}{4} \frac{d^2}{r^2}$$

$$\Omega_l = \Omega_s \iff r_s = \frac{d_s}{d_l} r_l$$
(2.1)

with diameters $d_l = 139$ mm for the large detector and $d_s = 56$ mm for the small one. The high-voltage bias for the two detectors were set to $U_l = 1400 \text{ V}$ for D_l and $U_s = 1000$ V for D_s . The gains of the amplifiers to 35 for D_l and 69 for D_s .

The signal from each detector is split into two branches which will be called the 'energy branch' and the 'time branch' (see figure 3). The modules in these branches are briefly characterized in section 2.2. The idea behind the 'time branch' is to record any coincident signal with the corresponding deviation from exact coincidence, which is expected to be gaussian distributed. Since the expectation value of the time delay between coincident events is zero and the time-to-amplitude converter (TAC) cannot record negative delays, the D_s -signal is shifted by a time that is large compared to the expected deviations. The 'energy branch' is supposed to eliminate random coincidents by opening the gate of the multichannel analyzer (MCA) for the 'time branch' signals only when the same event is measured by both detectors.



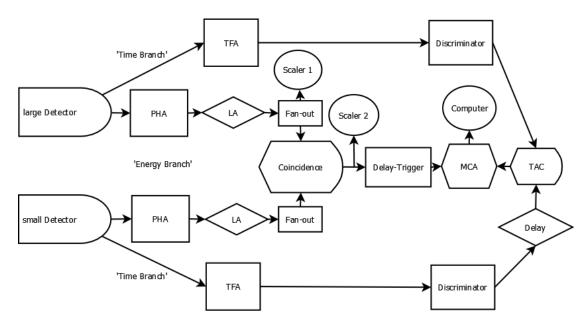


Figure 3: Flowchart showing all the electronic modules used for processing the detector signals

The 'time branch' starts at each detector with a timing filter amplifier (TFA) which is used to turn the detector signal into a clear analog signal. This signal is subsequently turned into a digital NIM pulse in a discriminator. The signal originating from D_l is forwarded into the start-channel of a TAC. The signal from D_s , on the other hand, is first delayed in a delay module by $t_0 = 48$ ns before being forwarded into the stop-channel of the TAC, which transfers a pulse to the ADC-channel of the MCA with an amplitude proportional to the time delay between the two incoming signals.

The 'energy branch' starts at each detector with an amplifier containing a pulse height analyzer (PHA). The signal produced by the detectors represents a broad energy spectrum (see section 3). Since only the γ -decays are of interest, the window of the PHA is set such that only the TTL pulses corresponding to the γ -ray emission are passed to the level adapter (LA), where they are converted into NIM signals. After the conversion, they pass a fan-out module where the signal of the D_l is split into two identical pulses and one of them enters a scaler S1 that counts the rate of events. This data is used for a cross-check of the measurement. The signals of the two detectors merge into an AND-gate inside a coincidence module. This module returns a pulse only if there is a signal coming from both detectors at the same time. The output signal is then forwarded into another scaler S2, which is also used for cross-checking, and into a delay-trigger that generates a rectangular TTL signal with adjustable delay and width. This signal is used as the gate input for the MCA.



The computer connected to the MCA records the number of events contained in a histogram where the bins correspond to the delay between the detector signals. This measurement is performed while placing the small detector at different angles θ relative to the large one, covering a θ -range from 180° to 60° with decrements of 15°. Each measurement will last for 20 min. In order to be able to assign an actual time delay to each bin in the histogram, a calibration measurement is carried out where the output of the TAC is observed for different known delays (see section 5.1).

2.2 Electronic Modules

A short description [1] for each electronic moudule used during the experiment is listed here, followed by two photographs showing the actual set-up, previously only illustrated schematically in figure 3.

The **pulse height analyser** (PHA) module creates a TTL output signal for every input signal inside a predifined range.

The **timing filter amplifier** (TFA) allows to adjust the shape of a signal according to the desired time constants.

The **discriminator** turns an analog input signal above a certain threshold into a logic (NIM) pulse whose width is adjustable.

The **level adapter** (LA) converts logic TTL signals into NIM signals or the other way around.

The fan-out module splits an input signal into several identical output signals.

The **coincidence module** contains a logic AND-gate that generates a NIM pulse whenever there is a signal on all input channels.

The **scaler** counts the number of incoming signals.

The **delay module** is a passive module that delays the incoming signal by a desired amount of time.

The **delay trigger** turns a rectangular input signal into an output signal with adjustable length and delay.

The time-to-amplitude-converter (TAC) returns a pulse whose amplitude is proportional to the time delay between the signal at the start channel and the signal at the stop channel

The **multichannel analyser** (MCA) takes pulses of different amplitudes as input and returns a histogram that counts the number of pulses per amplitude. It also contains a Gate channel which allows to register an event signal only when there is a signal on the Gate channel.

In this experiment, two different kinds of logic pulses are used. A TTL signal has a baseline of \pm 0.5V and is a rectangular pulse corresponds to a voltage of 5V. A NIM signal, on the other hand, is interpreted as a pulse when it falls to \pm 800mV.





Figure 4: Photograph of the actual modules during the experiment: 1)PHA of D_l , 2)TFA of D_l , 3)PHA of D_s , 4)TFA of D_s , 5)LA, 6)Fan-out, 7)Coincidence module, 8)MCA



Figure 5: Photograph of the actual modules during the experiment: 9)Fan-out, 10)Discriminator, 11)Delay module, 12)TAC, 13)Scaler S1, 14)Scaler S2, 15)Delay Trigger



2.3 Scintillation Counters and Photomultipliers

The detectors used in this experiment are scintillation counters. They consist mainly of a scintillating material that is optically coupled to a photomultiplier tube.

The molecules of a scintillating material (in this case NaI crystals doped with Thallium) become excited when radiation passes through and emit light during deexcitation (here with a wavelength of 410 nm). The intensity of this light emission depends linearly on the energy deposited. The light is then detected by a photomultiplier tube. [3]

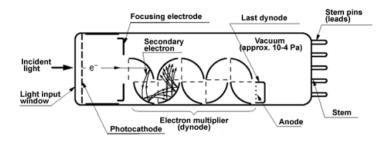


Figure 6: Schematic illustration of a photomultiplier tube [3]

The concept of a photomultiplier tube is illustrated in figure 6. When light above a threshold energy hits the photocathode, an electron is ejected from this photosensitive material due to the photoelectric effect. The electron is accelerated towards an array of electrodes calles dynodes where the potential is increasing step by step at every dynode. When an electron strikes a dynode, more electrons are ejected from the material which are then accelerated towards the next dynode. Thus, the amount of electrons increases exponentially until they hit the anode in the end of the tube where they induce a measurable current that can be used as an analog signal. The cathode and dynode systems are assumed to be linear, meaning that the current signal at the output is directly proportional to the incoming amount of photons. [3]



3 Energy Spectrum

Figure 7 shows the energy spectra measured at the two detectors correspondingly where the horizontal axis has been rescaled. The original scale can be found in figure 13 in the appendix, where the number of events was merely a function of ADC-channels. The two sharp peaks in the right half of the plot could be easily identified as the typical γ -decay absorption lines. Their corresponding energy values of $E_{\gamma,1}=1172$ keV and $E_{\gamma}=1332$ keV are well known in the literature and were used for the energy calibration that made the rescaling to energy units possible.

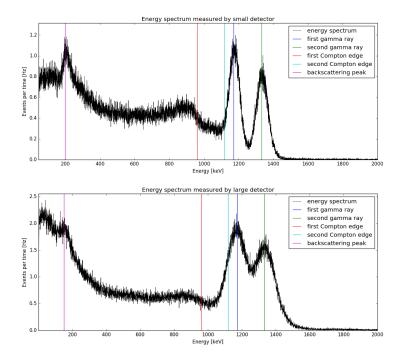


Figure 7: Energy spectra of each detector with tagged γ -ray absorption lines and typical detection byproducts

In addition to the two sharp gamma ray peaks, there is a typical Compton continuum visible on the left. The Compton edge can be obtained as follows [4]

$$E_{CE} = E_{\gamma} \cdot \left(1 - \frac{1}{1 + 2E_{\gamma}/m_e c^2}\right).$$
 (3.1)

This yields a value of $E_{CE,1} = 963$ keV for the first γ -ray and $E_{CE,2} = 1118$ keV for the second γ -ray. The first Compton edge is clearly visible in the plots while the second one coincides with the finite width of the first γ -ray line. Thus, it is clear that the Compton continuum does not vanish completely after the first edge.



There is also an additional peak visible in both plots at around 200 keV in the case of the small detector and 150 keV for the large one. This is a very common effect that can be traced back to backscattering which originates from γ -rays that have Compton scattered in the material surrounding the detector material before being detected. In the literature, this peak is said to be expected at around 200 keV, which is in very good agreement with the measurement. It also makes sense that this peak is not the same for both detectors, since they are not identical models.



4 Measurement Log

For each angle in the range of 60° to 180° with steps of 15° we took three 20 minute measurements: The number of coincident events measured with the software Maestro r_m , the number of coincident events measured by one of the scalers r_s , and the total number of events measured by the fixed detectors r_{tot} . The results of these measurements for each angle are listed in table 1.

Angle θ [°]	r_m [Hz]	r_s [Hz]	r_{tot} [Hz]
60	3.54	4.34	1470
75	3.35	4.09	1469
90	3.33	4.12	1469
105	3.39	4.16	1468
120	3.44	4.20	1464
135	3.52	4.29	1464
150	3.60	4.35	1468
165	3.57	4.35	1466
180	3.67	4.31	1453

Table 1: Measurement results

The software Maestro gave us not only the number of coincident events as listed in the table above, but the number of events for 8191 ADC channels. As an example, the histogram for $\theta = 90^{\circ}$ is shown in figure 8, the rest can be found in the appendix on page 21.

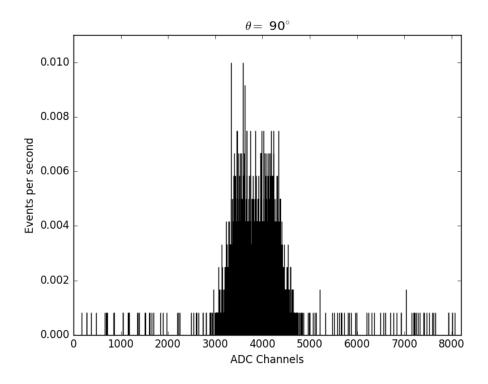


Figure 8: Histogram with ADC channels for $\theta=90^\circ$



5 Data Analysis

5.1 Time Calibration

Equivalently to the energy calibration in section [3], the MCA transmitted the number of events as a function of the corresponding channel to the computer. In order to convert the channel number to the corresponding physical time delay, a time calibration measurement has been performed. This consisted of splitting the signal of one of the detectors and then connecting one copy directly to the start-channel of the TAC and the other one to the stop-channel after passing the delay module first. The signal that the TAC induces for various delay settings in the MCA could then be directly observed in the computer. Data was collected and saved for five different delay settings, which can be seen in figure 9.

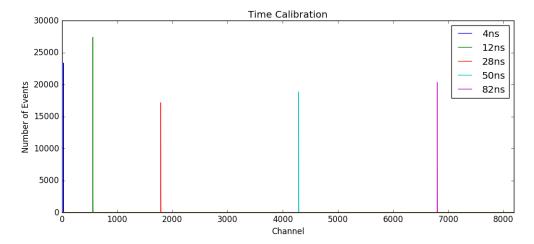


Figure 9: Time calibration measurement where each line corresponds to a known delay time

From this, we could deduce which channel corresponds to which magnitude of delay. The channel number and their corresponding delay time were plotted and a linear fit through the data points could be performed (see figure 10), which yields the following rescaling function where \hat{T} corresponds to the time-scale in [ns] and \hat{C} to the channel-scale

$$\hat{T} = \hat{C} \cdot 0.0111 \text{ns} + 5.2041 \text{ns} \tag{5.1}$$



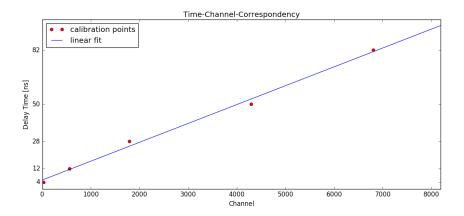


Figure 10: Linear fit through the corresponding delay-channel-points

5.2 Gaussian Fit

We started by rescaling the x-axis from ADC channels to time (see section 5.1). To make sure there are enough counts per bin we had to rebin the data. We merged 20 bins to one, giving us a time resolution of 0.22 ns. We then used the χ^2 -method to fit a Gaussian function on each histogram, taking into account an uncertainty on each bin of

$$\Delta_{bin} = \sqrt{n_{bin}} \tag{5.2}$$

Figure 11 shows the histogram with gaussian fit for the angle $\theta = 90^{\circ}$. The rest of the histograms can be found in the appendix on page 22.

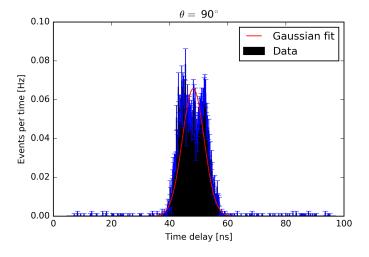


Figure 11: Histogram with gaussian fit for angle $\theta = 90^{\circ}$



We noticed that the measured data does not show a single peak at the set delay time, but instead two peaks slightly below and above. A reason for this might be the poor stability of the trigger settings in the discriminator. In a next step, we subtracted the background from the data. In order to estimate the background, we took the average of the bins in a range of $[3\sigma, 4\sigma]$ on both sides of the mean μ . Then we subtracted this average from each bin in the histogram. We then summed all the bins in a range of $[\mu-3\sigma, \mu+3\sigma]$ in the histogram to get the number of coincident events per second r_c for each angle. The statistical uncertainty was estimated by taking the square-root of the number of counts and normalizing it accordingly. Table 2 shows the result for each angle compared to the number of coincident events per second r_s as measured by the scaler.

Angle θ [°]	r_c [Hz]	r_s [Hz]
60	3.42 ± 0.05	4.34
75	3.24 ± 0.05	4.09
90	3.23 ± 0.05	4.12
105	3.30 ± 0.05	4.16
120	3.36 ± 0.05	4.20
135	3.42 ± 0.05	4.29
150	3.50 ± 0.05	4.35
165	3.47 ± 0.05	4.35
180	3.54 ± 0.05	4.31

Table 2: Rate for each angle

5.3 Angular Correlation Function

In order to be able to compare this data to the predicted value, we normalized the data to $r_c(\theta = 90^\circ)$. Then we used the χ^2 -method to fit a function in the form of equation 5.3 to the data points (figure 12).

$$W(\theta) = a_0 + a_1 \cos^2 \theta + a_2 \cos^4 \theta \tag{5.3}$$

This corresponds to formula 1.2 with leaving a_0 as a free parameter instead of constraining it to exactly 1.



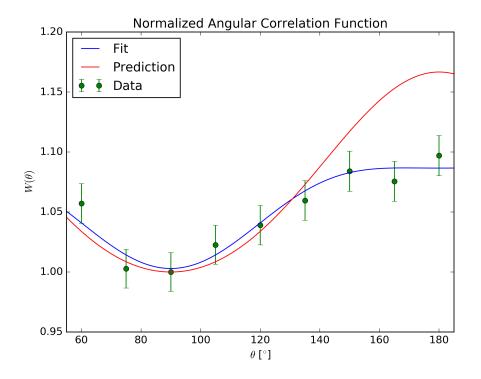


Figure 12: Normalized angular distribution

For the coefficients a_i we obtained:

$$a_0 = 1.003 \pm 0.007$$

$$a_1 = 0.175 \pm 0.044$$

$$a_2 = -0.091 \pm 0.043$$



6 Conclusion

6.1 Results

Comparing the results

$$a_0 = 1.003 \pm 0.007$$

 $a_1 = 0.175 \pm 0.044$
 $a_2 = -0.091 \pm 0.043$

derived in section 5.3 to the literature values [2]

$$a_{0,lit} = 1$$

 $a_{1,lit} = 1/8 = 0.125$
 $a_{2,lit} = 1/24 \approx 0.042$

we see that $a_{0,lit}$ is clearly within the range while $a_{1,lit}$ is slightly out of range. $a_{2,lit}$, on the other hand, is notably out of range. From figure 12 it is clearly visible that our measured data does not verify the prediction for angles starting at $\theta = 135^{\circ}$. In the range of $\theta \in [60^{\circ}, 120^{\circ}]$ the measured data corresponds to the prediction.

6.2 Discussion of Error Sources

A closer look at the vertical axes in figure 7 yields that the larger detector detector around twice as many events as the smaller one, despite the adjustments in behalf of solid angle coverage. For a better measurement, it might be wise to use identical detector models at equal distances.

The scaler S2 merely measures the amount of approximately coincident events from both detectors. Using the 'time-branch', more information is provided. Namely, the overlay of the TAC-signal and the gate-input from the delay trigger add more nuance by providing a quantitative delay between allegedly coincident events. As expected, the number of events counted by S2 was significantly higher at every angel. Nevertheless, it supplied a cross-check in the sense that the sclaer values behave the same way as the ones obtained by the TAC-signal. This similarity verifies that the electronic modules act as expected.

A significant source of uncertainty might come from the poor stability of the trigger settings in the discriminator that let to the splitting of the expected gaussian distribution into two discrete peaks. In order to avoid such error, one would have to be more careful



when adjusting the trigger in the future. Another way to deal with this problem would have been using a NIM signal again for the input of the discriminator.

In general, the measurement lacked a good degree of stability. A test run yielded that even slightly touching the cables could have an effect on the measurement. Not to mention the effect that the presence of the students and assistants might have had on the outcome. Thus we conclude that a more isolated setting could further improve the measurement.



7 Appendix

7.1 Energy Spectrum

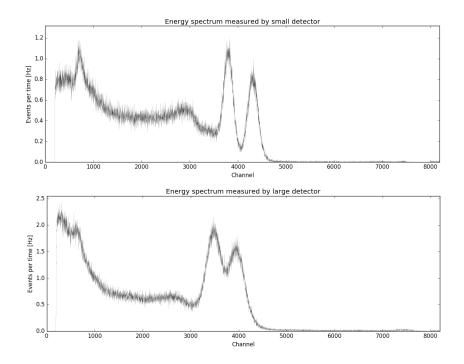


Figure 13: Energy spectrum measured as a function of ADC channels



7.2 Measurement Log

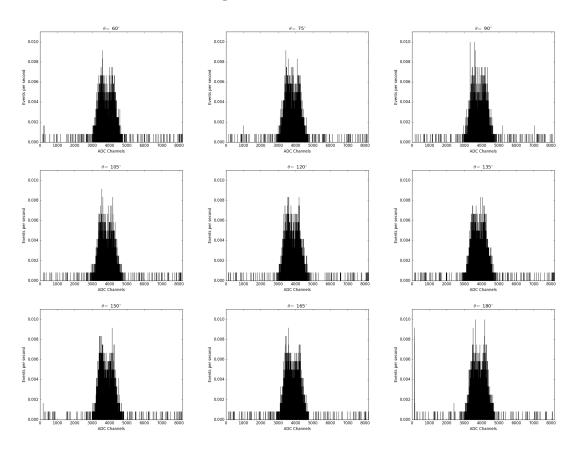


Figure 14: Histograms with ADC channels for all angles



7.3 Data Analysis

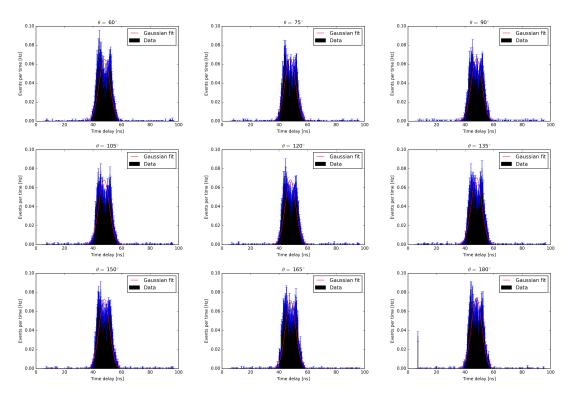


Figure 15: Rebinned histograms on time axis with background subtracted and fitted Gaussian for all angles

7.4 References

- [1]: Instruction manual and assistants
- [2]: Donald R. Hamilton Phys. Rev. 58, 122 (1940)
- [3]: Dr. William R. Leo Techniques for Nuclear and Particle Physics Experiments: A How-to Approach (1994)
- [4]: Glenn F. Knoll Radiation Detection and Measurement (2010)